

AIR TECHNICAL SETUP SUMMARY

[Last Updated: Aug 2025]

Objective: This user guide provides a summarized procedure for setting up AIR at a facility.

Prerequisite:

- Internet connectivity
- Accurate data / report


Linking up TaifaCare with OpenHIM

Steps:

- Log into TaifaCare as Admin



- Navigate to the developer tab. Then to legacy admin UI.

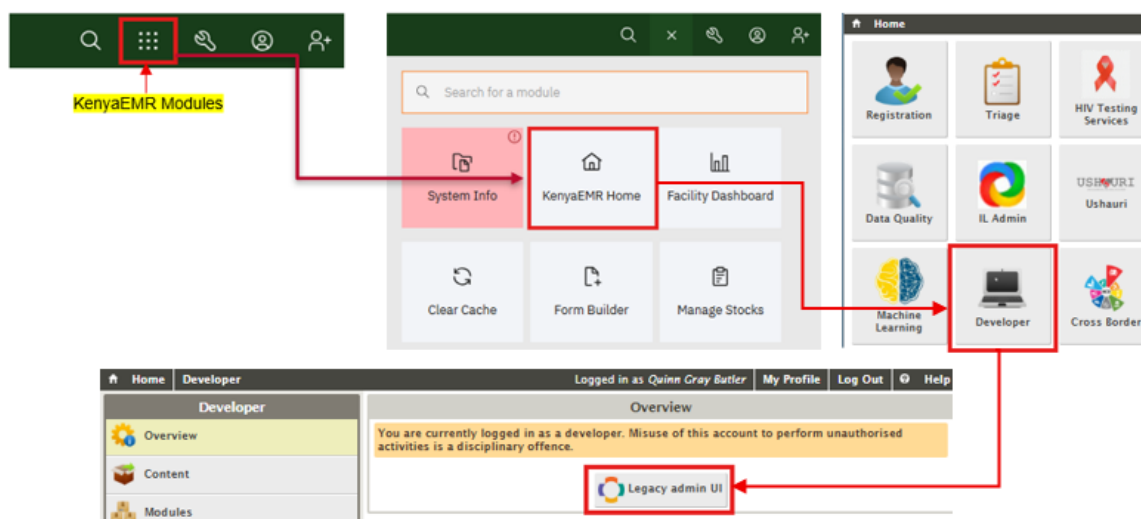

- Look for the maintenance part where you will see advanced settings.

Maintenance

[Set Implementation Id](#)

[System Information](#)

[View Quick Reports](#)

[Settings](#)

[Advanced Settings](#)

[View Server Log](#)

[View Database Changes](#)

- Once in the advanced settings search and fill the following global parameters:
 - o dhis.password and dhis.username (**Enter the correct credentials for a user with KHIS data entry rights for the respective site**)

dhis.password	kogutu	Remove
Password for DHIS server		
dhis.username	Kogutu@123	Remove
Username for DHIS server		

- **ilServer.address** : this is where OpenHIM outbox message will be routed

ilServer.address	https://openhimapi.kenyahmis.org/rest/api/IL/MOH_731/test	Remove
This is where IL Outbox messages will be routed to		

- **facility.mflcode** : facility MFL code. If it is not there then Scroll to the bottom of the page and click on add property.
- Two rectangular text boxes will appear.
- On the text box that is on your left type in **facility.mflcode** and on the right enter the actual MFL code.

facility.mflcode	13258
------------------	-------

Test the configurations:

- Log into TaifaCare
- Extract MoH 731 Report for the last month
- View report as ADX format
- Check to see if the IL Server address is prefilled as shown below.
- Try submitting data.

ADX Message for MOH 731

Reporting Date
Start Date: 2022-03-15 End Date: 2022-06-15
Server Settings
IP Address/URL: http://localhost:9721/api/
Edit

NOTE:

If you get error while submitting the data, this can be as a result of:

- Lack of internet connectivity
- Wrong KHIS login credentials
- Wrong DHIS end-point configurations
- Wrong headers on the ADX payload (this is very rare)

Check through the above as you eliminate each.

Further support

If in need of additional support, please contact us on the following toll-free line: **0800722440**

THE END